


DV-E™ Viscometer

our lowest cost digital viscometer

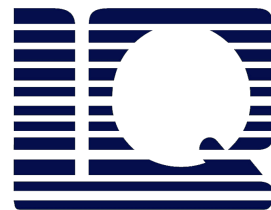

LABEQUIM, S.A DE C.V

No calculations required

- Direct reading of viscosity in cP or mPa-s

Displayed Info:

- Viscosity (cP or mPa-s)
- % Torque
- Speed/Spindle

Easy to Use:

- Flip a switch
- Turn a knob

Spindle/Speed Selection

Flip to "Speed"

- Turn the knob
- Choose the RPM

Flip to "Spindle"

- Turn the knob
- Choose the spindle

Auto Range

push for determining full scale range (FSR) viscosity

18 speeds for complete range capability

Accuracy: ±1.0% of range

Repeatability: ±0.2%



What's Included?

Instrument

6 spindles (RV/HA/HB) or 4 spindles (LV)

Spindle Guard Leg*

Lab Stand (Model A)

Carrying Case

*Not applicable to HA or HB models

Optional Accessories

Viscosity Standards

RV/HA/HB-1 Spindle

Quick Action Lab Stand

Temperature Bath

Small Sample Adapter

UL Adapter

Thermosel

Helipath Stand with T-bar Spindles

Spiral Adapter

DIN Adapter

Quick Connect/Extension Links

Vane Spindles

MODEL	VISCOSITY RANGE cP(mPa•s)		SPEEDS	
	Min.	Max.	RPM	Number of Increments
LVDV-E	1†	2M	.3-100	18
RVDV-E	100††	13M	.3-100	18
HADV-E	200††	26M	.3-100	18
HBDV-E	800††	104M	.3-100	18

† 1 cP achieved with UL Adapter accessory. 15 cP on LV with standard spindles.

†† Minimum viscosity is achieved with optional RV/HA/HB-1 spindle.
M = 1 million cP = Centipoise mPa•s = Millipascal•seconds



LABEQUIM, S.A. DE C.V.

labequim.com.mx

ventas@labequim.com.mx

TEL/FAX: 01 (222) 219 8815, 274 5072, 274 5073

25 ote. #217 San Pedro Cholula, Puebla, Pue; Col. San Pedro Mexicaltzingo